

# SOFRELINER TOUGH®

Soft Denture Reline Material

**Sofreliner Tough®** is an addition-cured silicone chairside soft lining material for dentures. Sofreliner Tough® has no odor, taste or chemical irritations. It offers a high tear resistance, high stain resistance and outstanding adhesion to the denture base. Sofreliner Tough® is designed to be used when durability is the primary concern.

The paste cures in 5 minutes intraorally, and the kit comes with 2 different finishing points (Shape adjusting point and Finishing point).




## Benefits

- Durability
- Strong adhesion
- High tear resistance
- High stain resistance
- High wear resistance
- Available:
  - Special burs for polishing and finishing
  - Silicone Remover (sold separately) for the removal of old reline materials
  - Step-by-step instructional video
- Auto-dispensing system/No hand mixing

## Indications

- Pain caused by sharp ridges of alveolar bone and atrophy of mucosa
- Poor retention caused by extreme ridge atrophy
- Major undercuts or maxillary/mandibular tori
- Retention of over dentures before the final attachment is placed

### Standard Handling Time for Sofreliner Tough®

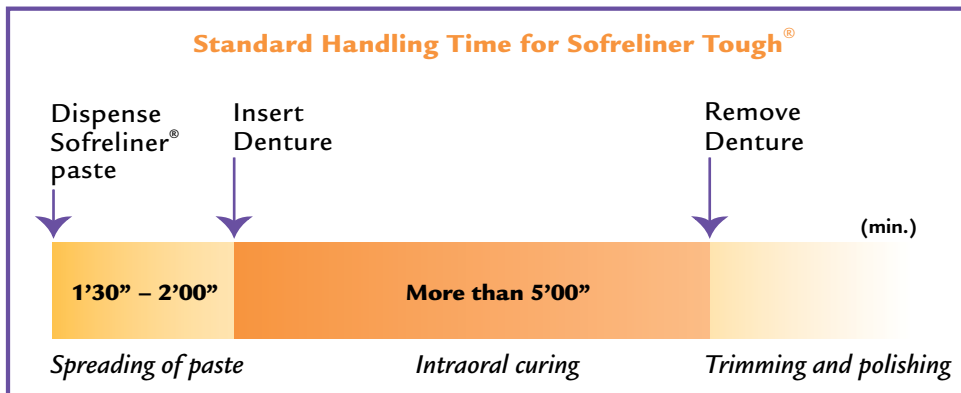

For more information, please call 1.877.ESTELITE or visit [www.tokuyama-us.com](http://www.tokuyama-us.com)



# SOFRELINER TOUGH®

Soft Denture Reline Material

## Product Packages Available

### Sofreliner Tough® Kit

- Paste/Base paste (27g), catalyst paste (27g)
- Primer (10ml)
- TD II Mixing Tips XS (x 10)
- Finishing point (Brown)
- Shape adjustment point (White)
- Diamond round bur
- Spatula
- Drip-cap for primer
- Plastic cup (x 2)
- Brush no. 5 (x 2)
- Marking pen



### Refills

- Sofreliner Tough® Paste/Base paste (27g), catalyst paste (27g)
- Sofreliner Primer (10ml)
- TD II Mixing Tips XS (50/box)
- TD II Mixing Tips S (50/box)
- Finishing point (Brown)
- Shape adjustment point (White)
- Diamond round bur
- Silicone Remover (7ml)
- TD II Dispenser



### Silicone Remover

The Silicone Remover (sold separately) is now available, which removes old silicone denture liner from the acrylic denture easily.



### SOFRELINER TOUGH® PRODUCTS

ITEM #	ITEM DESCRIPTION
23309	Sofreliner Tough® Kit
23326	Sofreliner Tough® Paste
23346	Sofreliner Tough® Primer
23366	Sofreliner Tough® Shape Adjustment Point
23367	Sofreliner Tough® Finishing Point
23368	Sofreliner Tough® Diamond Round Point
23337	Silicone Remover
23380	TD II Mixing Tip XS
23381	TD II Mixing Tip S
23372	TD II Dispenser

For more information, please call 1.877.ESTELITE or visit [www.tokuyama-us.com](http://www.tokuyama-us.com)

